

STRAIGHT STAIR LIFTS

A Smooth Step to Independence



Pinnacle



Pinnacle Heavy Duty



Vantage



Summit Outdoor



Harmar

Harmar is a leading manufacturer of lifts and is inspired by those who are mobility challenged. Harmar is driven to create mobility and accessibility products that help or enable people to go where they want, when they want, how they want.

Founded in 1998, Harmar has a proud history of growth by evolving our product lines to enrich the lives of our customers. We give people the path to independence to remain in their home, enjoy their families, and be engaged in their communities. Through innovation we are committed to help people of all ages overcome their limited physical mobility at home and away. All Harmar lifts are made in the USA, which means quality you can depend on.

Enjoy Your Home on All Levels

Whether due to an injury or illness, age or ability, the thought of getting up and down a staircase becomes more challenging. Safety and security are primary concerns, but so is being able to enjoy your life...your family. Don't let your staircase become an obstacle to your happiness.

Innovative engineering and design from a trusted manufacturer – Harmar, gives you the confidence and freedom to go wherever you want in your home. Whenever you want. All Harmar lifts are made in the USA, which means quality you can depend on.

Every year nearly 100,000 people invest in their independence by purchasing a stair lift. You can too. Be independent. Be confident. Be informed about your options.

Harmar is focused on developing products to help you. Our selection of straight stair lifts can handle the unique demands of a variety of homes.



90%

of people over Age 65 want to stay in their home for as long as possible. 19 million people with disabilities in the US alone have problems walking or climbing stairs.

The Pinnacle SL600

The Pinnacle SL600 straight stair lift is engineered to exceed the highest industry standards. Designed for optimal comfort, convenience and aesthetics. The most compact stair lift on the market—Pinnacle offers a space-saving 10.5” folded width when not in use. This lift has a capacity of 350 lbs and blends with any décor. Harmar’s patented nylon-polymer worm drive offers a smooth ride and requires no lubricants.

Smooth, Stable Ride

Compact Profile

Energy Efficient



Safety sensors on chassis and footrest instantly stop the lift at any obstacle.



Patented helical worm gear, requires no messy greases to maintain a smooth ride.



Two wireless remotes standard.



Narrowest lift on the market, fits in tight stairwells and makes it easy for others to pass.

Premium Stair Lift



Pinnacle Folding Rail

The **Pinnacle SL600FR Folding Rail** is an option that eliminates an obstruction at the bottom of the stairs. It is fully automatic yet 100% mechanical. This exclusive Harmar design option eliminates the need to stop and push a button or leave the seat.

- Great for doorway applications
- Valued option for busy families to eliminate trip hazards
- Protective rubber feet protect hardwood floors
- Fully automatic, no buttons or levers
- Charge strips above and below the hinge allow for continuous charge



Pinnacle Heavy Duty

The **Pinnacle SL600HD Heavy Duty** straight stair lift has a weight capacity of 600 lbs. High backed, contoured 25" wide seat ensures comfort and support. Sturdy, stationary armrests and reinforced footrest adds stability for safer transfers. It folds to 17" from wall when not in use. The Pinnacle Heavy Duty - unequalled lifting capacity.



Cost Effective Value of a Vantage SL400

Comfortable

Cost Effective

Confidence...



The Vantage SL400 straight stair lift is a great choice for those looking for a cost effective option. At 300 lbs capacity this sturdy, low maintenance and economical choice will provide a smooth ride ascending and descending. Complete with upholstered seat and powder coated chassis, this stair lift gets you up and down safely with style.

Peace of mind safety features include:

- Speed regulator maintains a comfortable, consistent speed
- DC drive system operates several cycles during power outages
- Footrest obstruction sensor, swivel seat, seat belt, and cut off switch



13.6" folded width for more room in your stairwell



Covered gear rack keeps your stairs clean




Two easy to use wireless remotes included

Summit Cable Drive Stair Lifts

The **Summit SL350AC** straight stair lift adjusts to your needs with three different settings for the backrest and footrest that can be adjusted to two heights for your utmost custom comfort.

- Operates on standard 115V AC household current, no batteries or charging required
- 350 lbs capacity
- Optional grocery basket

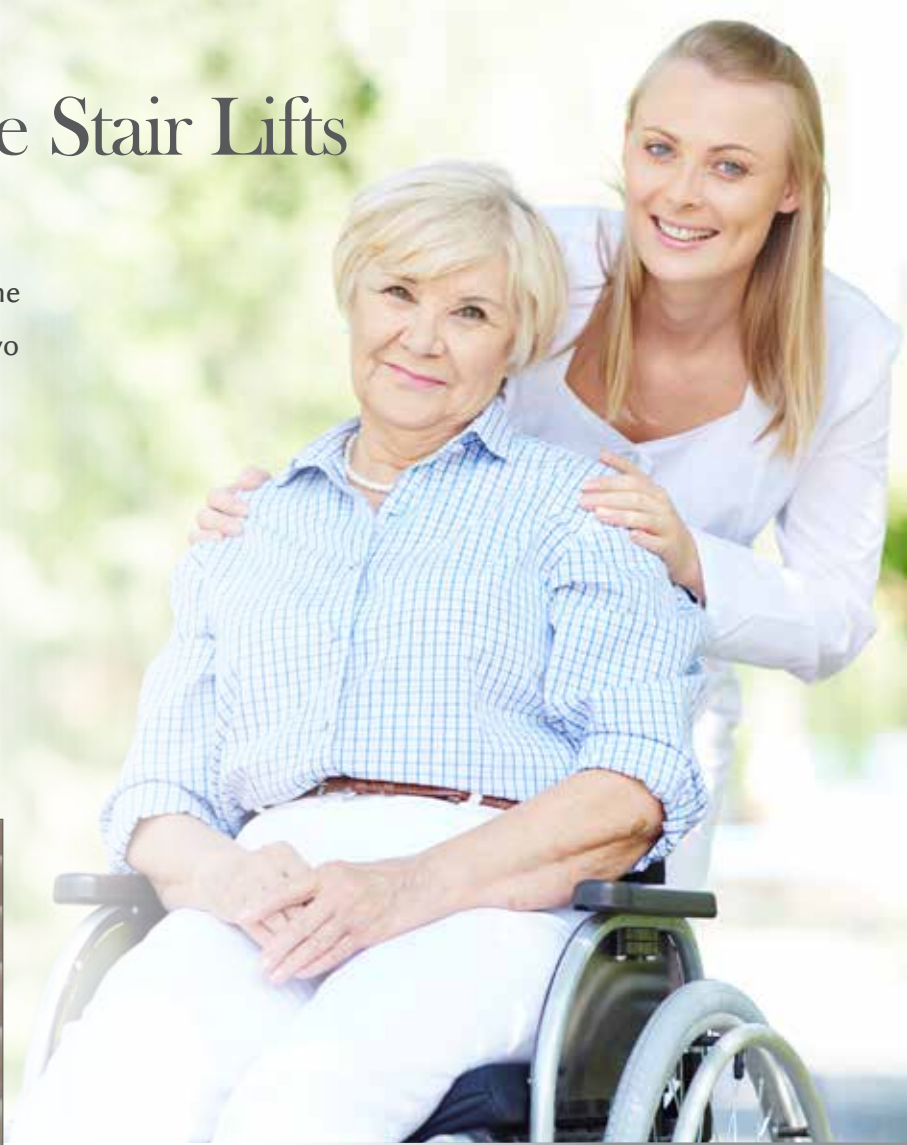

Summit SL350AC



14" width when folded



Wall mounted call/send controls



Safety sensors on footrest

Summit SL350OD for Outdoor Access

The **SL350OD Outdoor** straight stair lift is built and engineered to withstand the elements. The extruded aluminum track and marine-grade vinyl seating handle rain, wind, and snow.

- Designed for outdoor applications with operating temperatures ranging from 0° - 125°F.
- Zinc plated components
- Seat cushions in exterior marine grade vinyl
- Head cover, internal chassis cover, and large weatherproof cover included



SL350OD

SUMMIT

Harmar Straight Stair Lift Specifications

Model	Pinnacle SL600	Pinnacle SL600HD	Vantage SL400	Summit SL350AC	Summit SL350OD
Weight Capacity	350 lbs	600 lbs	300 lbs	350 lbs	350 lbs
Drive System	Nylon Polymer Helical Worm Gear	Nylon Polymer Helical Worm Gear	Steel Rack & Pinion	Cable Drive	Cable Drive
Speed	20 FPM	14 FPM	19 FPM	20 FPM	20 FPM
Seat	19"w x 16"h x 18"d	25"w x 16"h x 18"d	19"w x 16"h x 18"d	19"w x 14"h x 14"d	19"w x 14"h x 14"d
Folded Width	10.5"	17"	13.6"	14"	14"
Control Circuit	24V DC	24V DC	24V DC	24V AC	24V AC
Incline Limits	32° - 52°	32° - 52°	30° - 45°	30° - 45°	30° - 45°
Warranty	3-Year Parts 10-Year Gear Rack 1-Year Battery	3-Year Parts 10-Year Gear Rack 1-Year Battery	3-Year Parts 1-Year Battery	3-Year Parts	3-Year Parts
Options	Folding Rail Key Lock Additional Track Over 16'	Key Lock Additional Track Over 16'	Additional Track Over 16'	20' Track 25' Track (250 lbs) Key Lock 21" Wide Seat Basket	20' Track 25' Track (250 lbs) Key Lock 21" Wide Seat Basket

Harmar Offers Lifts for Home and Away



Curved Stair Lifts



Vertical Platform Lifts



Multiple auto lift solutions for wheelchair or scooter



Incline Platform Lifts



2075 47th St Sarasota, FL 34234
800.833.0478 www.harmar.com

*Images used for this brochure are for illustrative purposes only, please refer to local code for specific configuration requirements.